



Manual

Urinary Tract EKEGAS – OSR for BD MAX™
Version 03




400-028-C-MAX



24 reactions

For *In Vitro* Diagnostic Use
For use with BD MAX™ Open System Reagents on the BD
MAX™ System

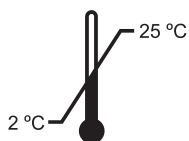

BioGXBV

Science Park 408, 1098 XH Amsterdam, The Netherlands
Phone: +31.20.893.4261
Fax: +31.20.240.9149

PROPRIETARY NAME

BioGX Urinary Tract EKEGAS – OSR for BD MAX™

INTENDED USE

The BioGX Urinary Tract EKEGAS – OSR for BD MAX™ is an automated *in vitro* diagnostic test reagent. The open system reagent (OSR) is used for the multiplex qualitative detection of DNA from *Escherichia coli*, *Klebsiella* spp., *Enterococcus faecalis*, Group A *Streptococcus* (EKEGAS), and a DNA sample processing control (SPC). The assay is run on an automated DNA extraction and real-time PCR instrument with urine from individuals at risk for presence of urinary tract infections. Automated extraction of the sample DNA is done using the BD MAX™ ExK™ DNA-3 series extraction kits. The extraction kits contain the SPC DNA so no external addition of SPC is required. The SPC serves as both an extraction control and an internal amplification control (IAC). Each tube of multiplex PCR mix is provided in BioGX proprietary Sample-Ready™ lyophilized format and contains all PCR primers, probes, enzyme, dNTPs, MgCl₂, buffers, and other components required for real-time PCR-based analysis of one sample.

SUMMARY AND EXPLANATION

Urinary tract infections (UTIs) are among the most common infections in people, and antibiotic treatment is usually helpful in treating an infection. UTIs can affect several parts of the urinary tract, but the most common type of UTI is a bladder infection (also known as cystitis). Infection can sometimes spread to the kidneys (known as pyelonephritis).

Women are at a higher risk for UTIs due to the proximity of their urethra to the anus as well as the having a shorter urethra than males. Most women will experience more than one UTI in their lifetimes. Other risk factors for UTIs include catheter use, previous UTIs, kidney stones, abnormal urinary tract (structure or function), and enlarged prostate.

Between 85% and 95% of uncomplicated UTIs are caused by infection with *E. coli*, however UTIs caused by a variety of other bacteria have been observed at lower rates including species from *Enterobacteriaceae*, *Enterococci*, *Klebsiella*, *Staphylococci*, Group A *Streptococci*, and Group B *Streptococci*.

PRINCIPLES OF THE PROCEDURE

The BioGX Urinary Tract EKEGAS – OSR for BD MAX™ is to be used with the BD MAX™ Open System for automated patient sample processing and molecular analysis. The BD MAX™ System uses a combination of lytic and extraction reagents to perform cell lysis and nucleic acid extraction. Following enzymatic cell lysis at elevated temperature, the released nucleic acids are captured by magnetic affinity beads. To control for extraction efficiency, a DNA Sample Processing Control is included in each BD MAX™ DNA Extraction Tube. The beads with bound nucleic acids are washed and the nucleic acids are eluted by heat in an elution buffer. The eluted nucleic acid is then mixed with the BioGX Rehydration Buffer, which is then transferred to the BioGX Sample-Ready™ lyophilized master mix tube in order to rehydrate the Sample-Ready™ lyophilized master mix. The rehydrated mix of amplification reagent and nucleic acid is then dispensed into the BD MAX™ PCR Cartridge. Microvalves in the BD MAX™ PCR Cartridge are sealed by the system prior to initiating PCR to prevent evaporation and amplicon contamination.

The amplified DNA targets are detected using hydrolysis probes labeled at one end with a fluorescent reporter dye (fluorophore) and at the other end with a quencher moiety. Probes labeled with different fluorophores are used to detect specific amplicons originating from *Escherichia coli*, *Klebsiella* spp., *Enterococcus faecalis*, Group A *Streptococcus*, and a Sample Processing Control in five (5) different optical channels of the BD MAX™ System: *Enterococcus faecalis* amplicons are detected in the 475/520 channel, Group A *Streptococcus* amplicons are detected in the 530/565 channel, *Klebsiella* spp. amplicons are detected in the 585/630 channel, *Escherichia coli* amplicons are detected in the 630/665 channel, and the DNA Sample Processing Control is detected in the 680/715 channel. When the probes are in their native state, the fluorescence of the fluorophore is quenched due to its proximity to the quencher. However, in the presence of their specific target cDNA, the probes hybridize to their complementary sequences and are hydrolyzed by the 5'-3' exonuclease activity of the DNA polymerase as it synthesizes the nascent strand along the DNA template. As a result, the fluorophores are separated from their quencher molecules and fluorescence is emitted. The amount of fluorescence detected in the five (5) optical channels used for the BioGX Urinary Tract EKEGAS – OSR for BD MAX™ is directly proportional to the quantity of the corresponding probe that is hydrolyzed, and therefore proportional to the amount of synthesized target. The BD MAX™ System measures these signals at the end of each amplification cycle in real time, and interprets the data to provide a qualitative result for each of the above targets.

REAGENTS

Qty	REF	Contents	Tests
2	400-028-MAX	BioGX Urinary Tract EKEGAS - OSR for BD MAX™ Sample-Ready™ lyophilized PCR Master Mix containing polymerase, nucleotides, specific molecular primers and probes, Sample Processing Control-specific molecular primers and probe.	12 tests per pouch
1	800-028-C	Rehydration Buffer Tube (C) Open System Reagents for BD MAX™ Reagent tube containing a rehydration buffer for use in Lyophilized PCR Master Mix rehydration.	24 tests per pouch

NOTE: Safety Data Sheets (SDS) are available at www.biogx.com/eu or by request.

EQUIPMENT AND MATERIALS REQUIRED BUT NOT PROVIDED

- BD MAX™ ExK™ DNA-3 (BD catalog no. 442822)
Extraction Kits include Sample Buffer Tubes (SBT), Septum Caps, Extraction Tubes, and Unitized Reagent Strips sufficient for 24 tests.
- BD MAX™ PCR Cartridges (BD catalog no. 437519)
- Urine Collection Device.
- Vortex Genie 2 Vortexer (VWR catalog no. 58815-234) or equivalent.
- Disposable nitrile gloves.

WARNINGS AND PRECAUTIONS



- Treat all biological specimens, including used Extraction Kits and PCR Cartridges, as if capable of transmitting infectious agents in accordance with safe laboratory procedures such as those described in CLSI Document M29¹ and in Biosafety in Microbiological and biomedical Laboratories.²
- Performance characteristics of this test have been established only with the specimen types listed in “Intended Use” section. The performance of this assay with other specimen types or samples has not been evaluated.
- Do not use the reagents if the protective pouches are open or torn upon arrival.
- Close reagent protective pouches promptly with the zip seal after each use. Remove any excess air in the pouches prior to sealing and store at 2-25 °C.
- Do not remove desiccant from the PCR Master Mix pouches.



- Do not use Master Mix if the desiccant is not present or is broken inside the Master Mix pouches.
- Do not use reagent tubes if the foil seal has been opened or damaged.
- Do not mix reagents from different pouches and/or kits and/or lots.
- Do not use expired reagents and/or materials.
- Each Master Mix and Rehydration Buffer tube is used to process a single sample. Do not reuse Master Mix or Rehydration Buffer tubes.
- Refer to BD MAX™ ExK™ DNA-3 Extraction Kit Instructions for information about proper handling, cautions, and proper waste disposal.
- Do not mix septum caps between Sample Buffer Tubes or re-use septum caps as contamination may occur and compromise test results.
- Check BD Unitized Reagent Strips for proper liquid fills (ensure that the liquids are at the bottom of the tubes).
- Do not pipette by mouth.
- Do not smoke, drink, or eat in areas where specimens or kits are being handled.
- Dispose of unused reagents and waste in accordance with country, federal, provincial, state, and local regulations.
- Use clean gloves when handling extraction kit components and PCR reagents and buffer tubes.

STORAGE AND STABILITY



- BioGX recommends long-term storage at 2-25°C.
- Reagents have been tested to demonstrate optimal performance when stored properly and consumed by the Expiration Date. Long-term stability studies are ongoing and the Expiration Date will be amended as additional data is available.
- Avoid exposing the reagents (lyophilized or rehydrated) to direct sunlight or long-term ambient lighting.
- Tightly reseal the pouch with unused reactions and immediately store the pouch in a dry location after opening.
- Avoid exposure to moisture and use the entire contents of the opened pouch within 1 month.

INSTRUCTIONS FOR USE

Install the BioGX Electronic User Defined Protocol on the BD MAX™

It will be necessary to import an Electronic User Defined Protocol (eUDP) onto the BD MAX™. The most current eUDP is available for download on www.biogx.com/eu by clicking on “Product Documentation” and selecting the appropriate platform and product name. eUDPs can also be obtained by emailing BioGX at eu@biogx.com. Please refer to the BD MAX™ user manual for uploading instructions.

NOTE: eUDPs are specific to extraction kit type and are programmed to be used with 3-snap strips unless otherwise indicated. If a 4-snap strip is used it is necessary to modify the eUDP to a 4-snap program by opening the eUDP in “Test Editor” and selecting the corresponding 4-snap extraction strip type in the “Extraction Type” drop down.

Specimen Collection/Transport

Urine specimens should be collected, transported, and stored according to institutional and laboratory standard operating procedures.

Specimen Preparation

Urine

Thoroughly vortex the sample prior to addition to the SBT. Add 500 µL of sample/specimen directly to the SBT.

Other Sample Types



This assay has been optimized for use with the sample types and volumes described above. Use of any other specimen type, collection method, or sample volumes may be inhibitory to the PCR or disrupt extraction without appropriate Guardrail and processing volume adjustments. BioGX does not make claims for processing methods or sample types other than those described in this product insert.

Setting up the Unitized Reagent Strip on the BD MAX™

1. Wear nitrile gloves when handling Sample-Ready™ lyophilized reagents to reduce the generation of static charges. **DO NOT** use latex gloves.



2. Use only BD MAX™ ExK™ DNA-3 extraction kits with the BioGX Urinary Tract EKEGAS – OSR for BD MAX™. DO NOT use BD MAX™ mastermix or the blank 0.3 mL conical tubes from the BD MAX™ ExK™ DNA-3 extraction kit.
3. Load one extraction cartridge into the extraction tray per specimen to be tested.
4. Snap one BD MAX™ ExK™ DNA-3 Extraction Tube into position 1 (Snap-in 1) of each Unitized Reagent Strip. (See Figure 1)

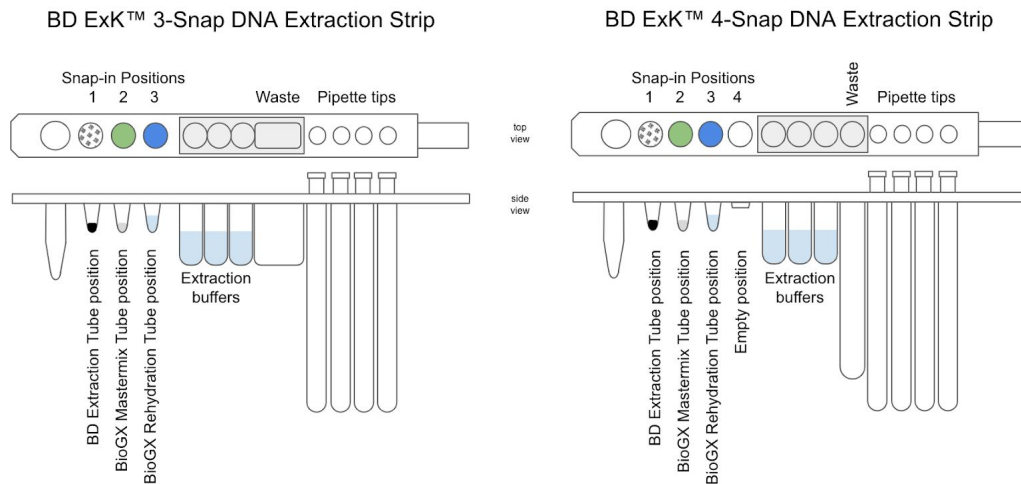


Figure 1 – Diagram of BD MAX™ ExK™ 3-snap and 4-snap Unitized Reagent Strips

5. Snap one BioGX Sample-Ready™ lyophilized PCR Master Mix reagent tube into position 2 (Snap-in 2) of each Unitized Reagent Strip. Check to make sure the Sample-Ready™ lyophilized cake is at bottom of tube prior to inserting into Unitized Reagent Strip. The funnel-shaped cake may be in any orientation (v, >, ^, <) in the **bottom** of the tube.
6. Snap one BioGX Rehydration Buffer tube into position 3 (Snap-in 3) of each Unitized Reagent Strip. Check to make sure the buffer is at bottom of tube prior to inserting into Unitized Reagent Strip.
7. Lift the tray and briefly examine the bottom of each Unitized Reagent Strip to ensure all reagents are at the bottom of each tube.
8. Proceed with worklist generation and sample loading per BD MAX™ operating instructions. Select the appropriate User Defined Protocol (eUDP) provided by BioGX. eUDPs are specific to extraction kit type and are programmed to be used with 3-snap strips unless otherwise indicated. If a 4-snap strip is used it is necessary to modify the eUDP to a 4-snap program by opening the eUDP in “Test Editor” and selecting the

corresponding 4-snap extraction strip type in the “Extraction Type” drop down.

9. Load the extraction tray and, if necessary, a new PCR card into the instrument, close the door, and click “Start Run.”

NOTE: Always first insert all Snap 1 tubes, then all Snap 2 tubes, then all Snap 3 tubes into the Unitized Reagent Strip.

NOTE: If using a 4-snap extraction strip, snap-in position 4 will remain empty.

QUALITY CONTROL

CONTROL

Each BioGX Urinary Tract EKEGAS - OSR for BD MAX™ includes molecular primers and probes specific for the detection of the DNA sample processing control (SPC) present in the BD MAX™ ExK™ DNA-3 Extraction Kit. No external addition of SPC is required. The SPC serves as both a sample extraction control and a PCR internal amplification control (IAC).

RESULTS INTERPRETATION

Results are available on the *Results* tab in the *Results* window on the BD MAX™ System monitor. The BD MAX™ System software automatically interprets the test result when the BioGX eUDP is used. Possible results for each target are shown in Table 1. Presence of one or more of the targets is possible, and will result in multiple targets being positive at once.

Results	Interpretation
Efaec POSITIVE	<ul style="list-style-type: none"> • The <i>Enterococcus faecalis</i> target has a Ct within the valid range and endpoint above the minimum setting.
GAS POSITIVE	<ul style="list-style-type: none"> • The Group A <i>Streptococcus</i> target has a Ct within the valid range and endpoint above the minimum setting.
KlebOP POSITIVE	<ul style="list-style-type: none"> • The <i>Klebsiella</i> spp. target has a Ct within the valid range and endpoint above the minimum setting.
Ecoli POSITIVE	<ul style="list-style-type: none"> • The <i>Escherichia coli</i> target has a Ct within the valid range and endpoint above the minimum setting.

Efaec NEGATIVE, GAS NEGATIVE, KlebOP NEGATIVE, OR Ecoli NEGATIVE	<ul style="list-style-type: none"> The respective target did not amplify and the SPC has a Ct within the valid range and endpoint above the minimum setting.
UNR	<ul style="list-style-type: none"> Unresolved Result. No target amplification; No SPC amplification.
IND	<ul style="list-style-type: none"> Indeterminate due to BD MAX™ System failure (with Warning or Error Codes*)
INC	<ul style="list-style-type: none"> Incomplete Run (with Warning or Error Codes*)

*Refer to the “Troubleshooting section of the BD MAX™ System User’s Manual for interpretation of warning and error codes.

NOTE: In the presence of a high concentration positive result for any target, the SPC may be adversely affected (no amplification or delayed). This is normal.

REPEAT TEST PROCEDURE

In case of instrument failure, repeat testing can be performed by setting up a new run using the original sample/specimen and a fresh SBT as described above in the Specimen Preparation section.

LIMITATIONS OF THE PROCEDURE

- This product is intended for use with specimens collected using specimen collection and transport devices listed in the “Equipment and Materials Required But Not Provided” section.
- This product should only be used with BD MAX™ Open System Reagents on the BD MAX™ System.
- Incorrect test results may occur from improper specimen collection, handling or storage, technical error, sample mix-up or because the number of organisms in the specimen is below the analytical sensitivity of the test. Careful compliance with the package insert instructions and the BD MAX™ System User’s Manual are necessary to avoid erroneous results.
- Good laboratory technique is essential for the proper performance of this assay. Due to the high analytical sensitivity of this test, extreme care should be taken to preserve the purity of all materials and reagents.
- A positive test result does not necessarily indicate the presence of viable infectious organisms. A positive result is indicative of the presence of

target nucleic acid. A negative test result does not preclude presence of infectious organisms and should not be used as the sole basis for treatment or other patient management decisions.

- As with all PCR-based *in vitro* diagnostic tests, extremely low levels of target below the limit of detection of the assay may be detected, but results may not be reproducible.
- False negative results may occur due to loss of nucleic acid from inadequate collection, transport or storage of specimens, or due to an inadequate cell lysis and/or extraction. The Sample Processing Control has been added to the test to aid in the identification of specimens that contain inhibitors to PCR amplification and as a control for reagent integrity and of the assay system as a whole. The Sample Processing Control does not indicate if nucleic acid has been lost due to inadequate collection, transport or storage of specimens, or if cells have been adequately lysed.
- The BioGX Urinary Tract EKEGAS – OSR for BD MAX™ results may sometimes be Unresolved due to an invalid Sample Processing Control, or be Indeterminate or Incomplete due to instrument failure, and require retesting that can lead to a delay obtaining final results.
- Mutations or polymorphisms in primer- or probe-binding regions may affect detection of new or unknown *Escherichia coli*, *Klebsiella* spp., *Enterococcus faecalis*, or Group A *Streptococcus* resulting in a false negative result with the BioGX Urinary Tract EKEGAS – OSR for BD MAX™
- The BioGX Urinary Tract EKEGAS – OSR for BD MAX™ requires the use of five (5) optical channels from the BD MAX™ System: 475/520 channel, 530/565 channel, 585/630 channel, 630/665 channel, and 680/715 channel.

PERFORMANCE CHARACTERISTICS

Analytical and Diagnostic Specificity

Specificity was determined by running negative sample matrix (Urine) spiked with positive control template. The BioGX Urinary Tract EKEGAS - OSR for BD MAX™ was positive for *Escherichia coli*, *Klebsiella* spp., *Enterococcus faecalis*, and Group A *Streptococcus*.

The BioGX Urinary Tract EKEGAS - OSR for BD MAX™ test was run with ATCC MSA-1002 (20 Strain Even Mix Genomic Material) which does not contain genomic DNA for *Klebsiella* spp. and Group A *Streptococcus* and does contain genomic DNA for *Escherichia coli* and *Enterococcus faecalis*. Results were negative for *Klebsiella* spp. and Group A *Streptococcus* and positive for *Escherichia coli* and *Enterococcus faecalis*.

The BioGX Urinary Tract EKEGAS - OSR for BD MAX™ test was run against the ZepTometrix NATtrol Blood Culture Panel (Product Reference: NATBCP-BIO) which contains intact bacteria from relevant species. The test successfully detected controls for *Enterococcus faecalis*, *Escherichia coli*, *Klebsiella oxytoca*, *Klebsiella pneumoniae*, and *Streptococcus pyogenes*. The test reported negative results for the following: *Acinetobacter baumannii*, *Candida albicans*, *C. glabrata*, *C. krusei*, *C. parapsilosis*, *C. tropicalis*, *Enterobacter cloacae*, *Haemophilus influenzae*, *Listeria monocytogenes*, *Neisseria meningitidis*, *Pseudomonas aeruginosa*, *Serratia marcesens*, *Staphylococcus aureus*, *S. epidermis*, *S. agalactiae*, and *S. pneumoniae*.

Analytical and Diagnostic Sensitivity

The analytical sensitivity for the BioGX Urinary Tract EKEGAS – OSR for BD MAX™ was determined as follows: Dilution series of positive synthetic DNA samples for each target were added to the SBT in duplicate. Analytical sensitivity (Limit of Detection, LoD) was defined as the lowest concentration at which 95% of all replicates tested positive.

Target	LoD (copies per SBT)
<i>Escherichia coli</i>	115
<i>Klebsiella</i> spp.	68
<i>Enterococcus faecalis</i>	76
Group A <i>Streptococcus</i>	163

Analytical sensitivity during co-infection was tested by challenging the BioGX Urinary Tract EKEGAS – OSR for BD MAX™ in pairs of high concentration (10,000X LOD) of one target against low concentration (1X LOD) of another for all possible pairs in the test. All low concentration targets were positive and were not outcompeted by amplification of the high concentration target.

Reproducibility

The reproducibility study was performed on *Escherichia coli* synthetic target template by three separate technicians independently on two BD MAX™ instruments. All users obtained equivalent results.

Manufacturing Reproducibility

Two independent lots were manufactured at and were found to be equivalent based on internally established QC acceptance procedures.

REFERENCES
















1. Clinical and Laboratory Standards Institute. Protection of laboratory workers from occupationally acquired infections; Approved Guideline. Document M29 (Refer to the latest edition).

2. Centers for Disease Control and Prevention and National Institutes of Health. Biosafety in microbiological and biomedical laboratories. Choosewood L.C. and Wislon D.E. (eds) (2009). HHS Publication No. (CDC) 21-1112.
3. Centers for Disease Control Website, www.cdc.gov, accessed April 9, 2018.
4. Mayo Clinic Website, www.mayoclinic.org, accessed April 9, 2018.
5. Wiedemann B, Heisig, A, Heisig P, Uncomplicated Urinary Tract Infections and Antibiotic Resistance - Epidemiological and Mechanistic Aspects.
6. Nicolle L et. al., Uncomplicated urinary tract infection in women. *Canadian Family Physician*, 2016, vol. 52 (pg. 612-618).
7. Edwards MS, Baker CJ. Group B streptococcal infections in elderly adults, *Clin Infect Dis*, 2005, vol. 41 (pg. 839-847)
8. Palocaren T, et. al. Beta haemolytic group A streptococci causing urinary tract infection. *British J Urology* Oct;74(4):444-6.

REVISION HISTORY

Version	Date	Description of Change
03	01 FEB 2019	Updated storage recommendations from 2-8°C to 2-25°C.
02	09 NOV 2018	Added use of BD ExK 4-snap
01	30 AUG 2018	Initial Release

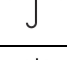SYMBOLS

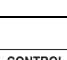Symbol	Meaning
	Catalog number
	<i>In vitro</i> diagnostic medical device
	Do not reuse
	Batch code
	Caution
	Consult instructions for use
	Manufacturer
	Contains sufficient for <n> tests
	Authorized Representative in the European Community
	Temperature limitation
	Keep dry
	Keep away from sunlight
	Expiration date
	Biological Risks
	Control



BioGXBV

Science Park 408, 1098 XH Amsterdam, The Netherlands
Phone: +31.20.893.4261 Fax: +31.20.240.9149